
Installing Log-Box Driver

1) Check whether the PC being used has 32-bit or 64-bit Windows installed.

Windows 8.1 Pro		
© 2013 Microsoft Corpor	ition. All rights reserved.	ows [.] 8
Get more features with a r		
stem		
Manufacturer:	Microsoft Corporation	
Processor:	Intel(R) Core(TM) IS-4300U CPU @ 1.90GHz 2.50 GHz	Surface
	4.00 GB	Surface
Installed memory (RAM):		

2) Go to <u>http://www.ftdichip.com/Drivers/VCP.htm</u> and download the appropriate driver. The example above would require the 64-bit driver.

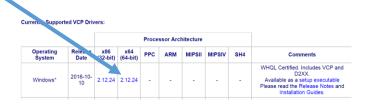

3) Unzip the driver to a suitable location.

v2.12.24 WHQL Certified	Y 🖉 Search CD	M v2.12.24 WHQL C F
Name	Date modified	Type
🗼 amd64	09/11/2016 10:13	File folder
👗 i386	09/11/2016 10:13	File folder
L Static	09/11/2016 10:13	File folder
ftd2xx.h	09/11/2016 10:13	H File
ftdibus	09/11/2016 10:13	Security Catalogue
(i) ftdibus	09/11/2016 10:13	Setup Information
ftdiport	09/11/2016 10:13	Security Catalogue
 ftdiport 	09/11/2016 10:13	Setup Information


4) Plug the USB to Serial cable (supplied with Logbox) into the PC. The cable does not need to be connected to the Logbox. The cable will be displayed in Device Manager (below).



5) Right click on Log-Box Cable FSN Logger and the select Update Driver Software.

Other devices	
Log-Box Cable FSN Log-	Update Driver Software
Standard Serial over Blu Standard Serial over Blu	Disable
Print queues Processors	Scan for hardware changes
Security devices	Properties

6) From the dialogue box that appears select Browse my computer for driver.

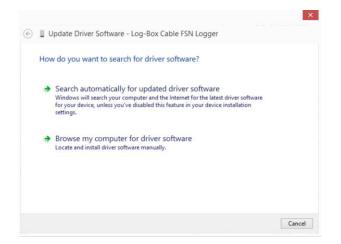

7) Click on the Browse button and navigate to the folder where the driver was unzipped in Step 3 above

0	te Driver Software - Log-Box Cable F! for driver software on your compute			
Search for	driver software in this location:		V Browse	
✓ Include	subfolders		 Browse 	
	Browse For Folder Select the folder that contains drivers for your ha	× rdware.		
⇒ L т я	Desktop Documents Dovuments Dovuments	^	mputer evice and all driver	
	CDM v2.12.24 WHQL Certified amd64 i386 i386		Next	Cancel
	Folder: CDM v2.12.24 WHQL Certified	>		
	ОК	Cancel	soft)	~

8) Click OK. Windows should then install the driver for the cable.

۲	Update Driver Software - USB Serial Converter	
	Windows has successfully updated your driver software	
	Windows has finished installing the driver software for this device:	
	USB Serial Converter	
	-	
		Close

9) Device Manager may show a USB Serial Port with an error (below)

Other devices
🕞 USB Serial Port
Ports (COM & LPT)
Standard Serial over Bluetooth link (COM3)
The standard Serial over Bluetooth link (COM4)

10) If the yellow ! triangle is shown right click on USB Serial Port entry and select Search automatically for updated driver software. Windows should update the driver automatically. A USB Serial Port will then be displayed in the Ports (COM & LPT) section of Device Manager.

4	Ports (COM & LPT)
	The standard Serial over Bluetooth link (COM3)
	🐺 Standard Serial over Bluetooth link (COM4)
	Transformation (COM6)
	-

The Log-Box can then be connected to the cable.

Installing Logbook Graphing

Logbook Graphing was updated in January 2016. The latest version is included with Log-box. If Log-Box was purchased before January 2016 the latest version can be downloaded from:


http://logbookgraphing.co.uk/product-category/tts/software/

Logbook Graphing requires the 32-bit version of Java. Please note: the 32-bit version of Java is required even if the PC is running 64-bit Windows.

If Java isn't installed when the Logbook is installed the error below will be shown.



Java can be downloaded from https://java.com/en/download/manual.jsp

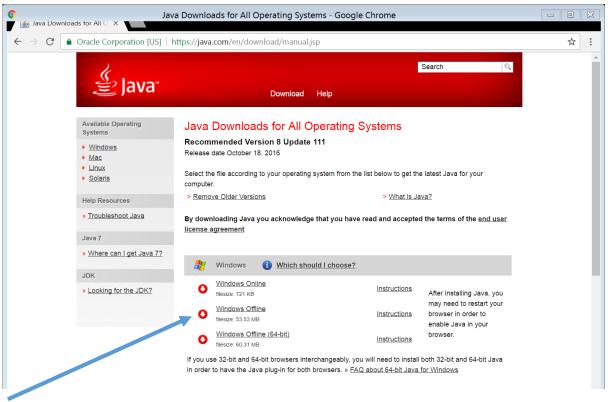

32-bit Java **must** be installed **before** Logbook Graphing is installed. If necessary uninstall Logbook, install 32-bit Java and re-install Logbook Graphing.